


NOTES

1. LED CORDSET PLUGS TO BE INSIDE CHAMBER ENVIRONMENT
2. PASS WIRES THROUGH MEMBRANE PLUGS IN RACEWAY COVER
3. LOOP CORDSETS INTO A KNOT INSIDE WIRE RACEWAY FOR STRAIN RELIEF

AMPERAGE DRAW: 0.3A @ 120VAC

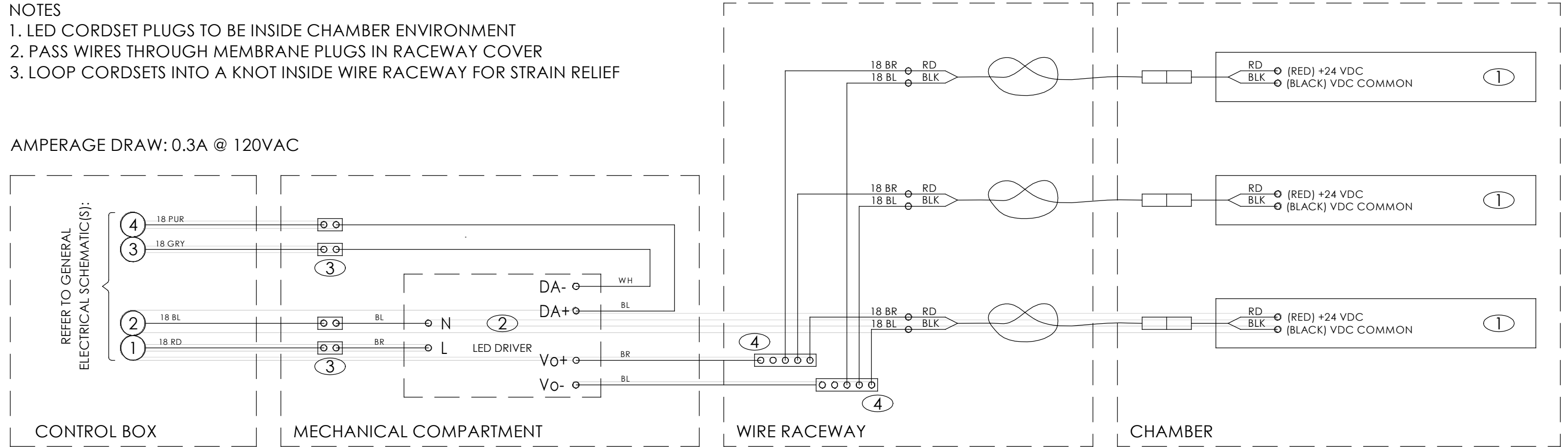

BOM:

- 1 - LED MODULE - 24" LENGTH, DIMMABLE, WHITE, 24 VDC, 5700K (ELD-106)
- 2 - LED DRIVER - 100-240, 277V INPUT, 24VDC OUTPUT, 50-60Hz, 120W (ELD-117)
- 3 - LEVER CONNECTOR - 2 TERMINAL, 24-12 GA WIRE, 600V, 20A (ECO-119)
- 4 - LEVER CONNECTOR - 5 TERMINAL, 24-12 GA WIRE, 600V, 20A (ECO-121)

TITLE: ELECTRICAL SCHEMATIC - GENERAL LIGHTING			
MODEL: ARC-30L-I	DRAWING:		REV:
BY: RAS	DATE: 6/5/23	032655-0-2	-